


seca 217

Stadiometer for mobile height measurement

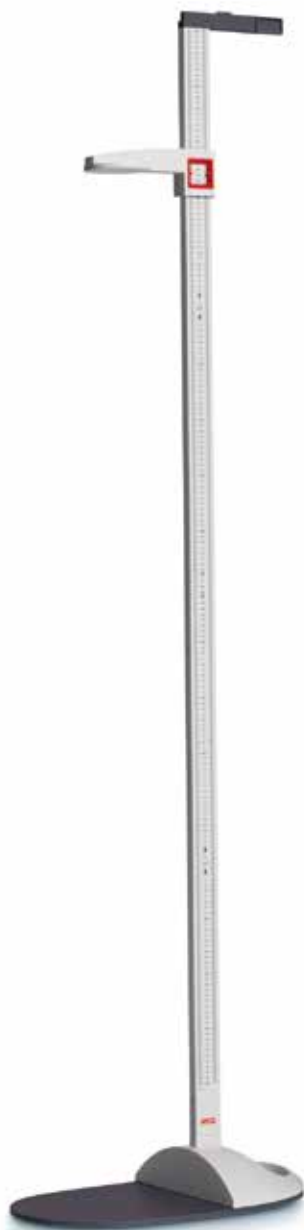

- Parts fit together perfectly for transport ease.
- Spacer for additional stability.
- Robust platform for secure stance.
- Adapter element seca 437 for connection to a flat scale.

seca 217:

Ideal solidity and stability.

The mobile stadiometer seca 217 is distinguished by its especially top-quality materials and a cleverly designed assembly system. The stability of this stadiometer is unmatched. It is particularly suited for use in doctors' practices and hospitals and is easily carried to schools for medical examinations and to patients' homes. When disassembled, the stadiometer can be conveniently transported to any location.



Fast and simple assembly without wall mounting.


The collapsible height measuring rod of the stadiometer seca 217 can be put together quickly and simply in just a few steps and attached to the stable platform. There's no need to attach the rod to the wall. It's easy to step onto the large base plate, which has no raised sides or edges.

Large base plate permits use anywhere.

The mobile stadiometer seca 217 can be used anywhere, thanks to its sturdy and large base plate. All the connections are very stable and the spacer between wall and rod prevent the rod from wobbling and trembling, movements which impair measurement accuracy. All these features make this stadiometer the ideal tool for anyone who needs a reliable measuring device for use at different locations.

Precise reading of results.

The large and stable head piece ensures highly accurate results. It is made of a firm non-warping plastic that smoothly glides into place. Another well-designed detail is found in the two-sided scale on the rod which permits a readout of results during the measurement process and guarantees precise measuring up to a height of 205 cm.

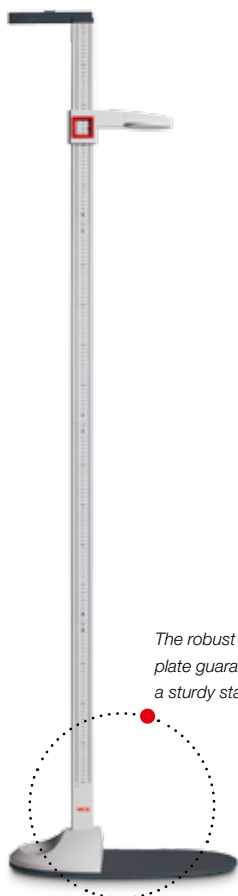

The robust base plate guarantees a sturdy stance.



The spacer provides extra stability and ensures precise measurement results.



The seca 437 adapter element connects the stadiometer to a flat scale.

To meet the challenges of daily medical practice.

The measuring rod can be taken apart in just a few steps and conveniently carried anywhere with the integrated handle. For stationary use the seca 437 adapter element joins the stadiometer to a flat scale, such as the seca 899, seca 878 or seca 877. The result is a weighing and measuring station that has proven itself in routine clinical work – as expected from the world market leader for medical weighing and measuring.

seca 217

Technical Data

- Measuring range: 20–205 cm, 20–205 cm / 8–81"
- Graduation: 1 mm, 1 mm / 1/8"
- Dimensions (WxHxD): 328 x 2,145 x 574 mm / 12.9 x 84.4 x 22.6"
- Weight: 3.6 kg / 7.9 lbs
- Optional: adapter element seca 437 (for combination with flat scale), carrying case seca 414 (for combination with flat scale and adapter element)

Exclusive Philippine Distributor:

WellnessPRO
tomorrow's medical devices

Unit 808 Richmond Plaza, San Miguel Ave cor Lourdes St.,
Ortigas Ctr. Pasig City 1605 Philippines
Mobile no.: +63 922 896 0964
Tel. no.: +63 2 636 3580, +63 2 637 0930, +63 2 637 0473
Fax: +63 2 636 7542
Email: wellnessproinc@gmail.com, jcjuan@wellnessproinc.com
www.wellnessproinc.com