

seca 417

Light and stable measuring board for mobile use



- Simple operation for quick mobile use.
- Folds up for space-saving storage.
- High-quality materials for long service life.
- Smooth surfaces for easy cleaning.

seca 417:

Heavy work made light.

The seca 417 measuring board is the mobile alternative to the seca 416 infantometer. Its fold-up mechanism and low weight make the seca 417 compact and lightweight for easy transport. It is simple to use and, thanks to its smooth surface, easy to clean. Durability and a long service life are guaranteed by high-quality materials.



Simple operation.

A hinge joins together the two parts of the seca 417 measuring board. It has an integrated head piece and a separate foot positioner that smoothly runs along the gently raised guide rail. For a precise measurement, all the user has to do is place the head of the child against the head piece, stretch out the legs and move the foot positioner to the soles of the feet.

Easy-to-read results.

The seca 417 stadiometer is ideal for mobile use. During regularly scheduled examinations, medical personnel check the nutritional condition and development of infants and toddlers. The process should go quickly, so it's important to be able to read results easily. The red markings and high-quality printing of the scale will not fade even after years of intense use.

Quick cleaning.

The smooth, rounded surfaces prevent bumps and jolts during measuring and make cleaning easy. The foot positioner can be removed to make the guide rails accessible for cleaning. The skin-friendly ABS plastic and the printed scale are impervious to all commercially available cleaning products and disinfectants.

Long service life.

The lightness of the measuring board is matched by the high quality of the folding mechanism, head piece and foot positioner. Special attention is paid to the choice of materials which can stand up to years of tough use without losing their shape. Additional protection is provided by the seca 412 carrier bag made of water-repellent lightweight nylon.




The measurement results are easy to read from the large scale.



The high-quality folding mechanism guarantees a long service life.

The sturdy head piece helps to position the child.

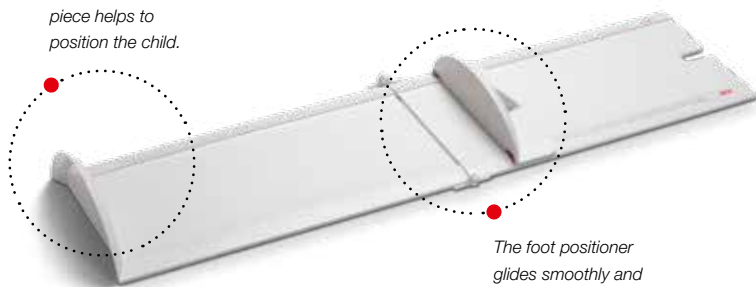

The foot positioner glides smoothly and precisely along special guide rails.

seca 417

Technical Data

- Measurement range: 10–100 cm, 10–100 cm / 4–39"
- Graduation: 1 mm, 1 mm / 1/16"
- Dimensions (WxHxD): 1,110 x 115 x 333 mm / 43.7 x 4.5 x 13.1"
- Dimensions when folded (WxHxD): 577 x 115 x 333 mm / 22.7 x 4.5 x 13.1"
- Device weight: 1.6 kg / 3.5 lbs
- Optional: carrying case seca 414 and seca 412

Exclusive Philippine Distributor:

WellnessPRO
tomorrow's medical devices

Unit 808 Richmond Plaza, San Miguel Ave cor Lourdes St,
Ortigas Ctr. Pasig City 1605 Philippines
Mobile no.: +63 922 896 0964
Tel. no.: +63 2 636 3580, +63 2 637 0930, +63 2 637 0473
Fax: +63 2 636 7542
Email: wellnessproinc@gmail.com, jcjuan@wellnessproinc.com
www.wellnessproinc.com