

seca 674

Electronic platform scales



- Generously sized platform and very high capacity.
- Suitable for weighing persons in a wheelchair or on a chair.
- Flexible cable remote display.
- Transport castors.

seca 674:

Weighing without restrictions.

Special situations demand special scales. Not every patient can step onto a scale quickly and easily. Most of all, patients with limited mobility or in wheelchairs and weak or seriously overweight persons need the help of a convenient scale like the seca 674. The platform scale serves a wide range of patients and eases the routine work of doctors and nurses.

Large platform with high load-bearing capacity.

Simply weigh the wheelchair too: The oversized non-slip platform of the seca 674 offers plenty of space for all conventional wheelchair models. For patients who have trouble standing, a chair can be placed on the scale – comfort and convenience that cannot be weighed. The Pre-TARE function simply subtracts the added weight of the chair or wheelchair and provides the weight of the patient alone. Thanks to its high load-bearing capacity of 360 kg, the scale is perfectly suited for use by seriously overweight persons. The ramp seca 470 is included in standard delivery.

Versatile wired remote display.


Because the display element of the seca 674 is connected to the scale by means of 2.7-meter-long cable, it can be positioned where the user can most easily read the results. A wall mounting is also included in the standard delivery. An additional feature is the RS232 adapter set seca 460 (optional) which connects this scale to a printer or PC.

Optional: non-slip base.

The seca 472 stand is the stable solution for cable remote displays of seca 674. The stable construction with non-slip rubber feet replaces the former model on transport castors and can be positioned everywhere quickly and securely. The display element can easily be attached to the stand by means of a stable bracket. The height of the display is continuously adjustable up to 100 cm.

Easy to transport.

Despite its size, the seca 674 is amazingly easy to handle. When upright, the scale can be moved easily on castors and steered by means of the integrated handle. This feature makes the seca 674 suitable for mobile use.

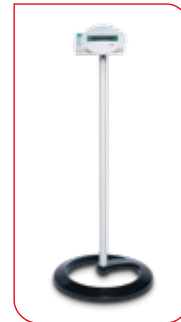

The seca 472 stand for cable remote displays of seca scales and measuring rods.



A ramp for easy roll-on of wheelchair is included in delivery.

seca 674

Technical Data

- Capacity: 360 kg / 800 lbs
- Graduation: 100 g / 0.2 lbs
- Dimensions (WxHxD): 904 x 70 x 1,064 mm / 35.6 x 2.8 x 41.9"
- Dimensions/platform (WxHxD): 760 x 55 x 930 mm / 29.9 x 2.2 x 36.6"
- Weight: 27 kg / 59.5 lbs
- Power supply: power adapter
- Functions: TARE, pre-TARE, auto-HOLD, BMI, kg/lbs switch-over, damping
- Optional: stand for cable remote display seca 472, second ramp seca 470, RS232 interface seca 460

Exclusive Philippine Distributor:

WellnessPRO
tomorrow's medical devices

Unit 808 Richmond Plaza, San Miguel Ave cor Lourdes St.,
Ortigas Ctr. Pasig City 1605 Philippines
Mobile no.: +63 922 896 0964
Tel. no.: +63 2 636 3580, +63 2 637 0930, +63 2 637 0473
Fax: +63 2 636 7542
Email: wellnessproinc@gmail.com, jcjuan@wellnessproinc.com
www.wellnessproinc.com